


PRODUCT DATA

I-Product Description & Classification		
I/1	Designation	Wood applicator
I/2	General Description	Applicator with Wooden Stick, 150x2.2mm, small cotton head
I/3	GMDN N°	34954
I/4	Reference Number	1102246
I/5	Class	<input type="checkbox"/> I <input checked="" type="checkbox"/> Is <input type="checkbox"/> IIa According to Directive 93/42/EEC - Annex IX, Rules 4 and 5
I/6	Notify Body	CE ₁₆₃₉
I/7	Composition Raw Material	100 % cotton and wooden stick
I/8	Sterile	<input checked="" type="checkbox"/> Yes <input type="checkbox"/> No If, yes <input checked="" type="checkbox"/> ETO <input type="checkbox"/> Gamma <input type="checkbox"/> Steam
I/9	First Packaging	Number of pieces in Packaging: 2 piece per peel pack
I/10	Picture	

II- Description of use and operating of device		
II/1	Intended use/indication:	A device intended used to apply a substance (eg medication) to the Ear/Nose/throat (ENT) treatment and to care the superficial wound for absorption, barrier or compressing Alone, single use and transient Use (<60min)
II/2	User:	Medical staff or professional
II/3	Contraindication of use	Not use, if the device is dirty or visual defects If packaging is damaged or open do not use. Sterile unless package is damaged.
II/4	Special requirement for use	Should not be use with Alcohol
II/5	Handling and storage conditions:	Standard conditions or atmospheric conditions

III- Supplementary information		
III/1	Compatibility with others medical devices (in the set)	n/a
III/2	Device contents animal substances	NO
III/3	Device contents medicinal substances	NO
III/4	Device comes into contact with substances	Applimed has demonstrated compatibility with the following substances: - Water - Peroxide water (3%) - Betadine - Octenisept - Eosine - Merfen
III/5	Used with Accessories	NO
III/6	Used in combination with other Medical Device	NO

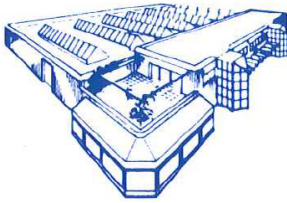


Zone industrielle
 Route de Pra de Plan 1
 CH-1618 CHATEL-ST-DENIS
 Switzerland

TEL : +41 (0) 21 948 92 74
 FAX : +41 (0) 21 948 90 36
 E-mail : info@applimed.ch
 Website : www.applimed.ch

IV- MANUFACTURING

IV/1	Manufacturing Operations	<input checked="" type="checkbox"/> peel pack <input checked="" type="checkbox"/> pouches <input type="checkbox"/> boxes <input checked="" type="checkbox"/> cartons <input type="checkbox"/> cleaning <input checked="" type="checkbox"/> sterilization	2 pieces per peel pack 100 peel packs per pouch 15 pouches per carton
IV/2	<p>Labelling: According to Directive 93/42 CEE Annex I § 13.3 and ISO 15223-1</p> <ul style="list-style-type: none"> Name and description of product Name and address of manufacturer LOT N° Article N° Year of manufacturing Expiry date Quantity Sterile Method of sterilisation Indication for single use Not use if the pack is opened CE 0120 	<p>Example of a label</p> <div style="border: 1px solid black; padding: 5px;"> <p>REF 1102246 LOT 21.04.27</p> <p>Code 11-2246 2026-03</p> <p style="text-align: center;">2x15x100 Stk/pieces</p> <p>Applicateur Bois 150x2.2mm Wattestäbchen steril / kleiner Kopf Cotton sticks sterile / small head Applikator sterile / lille hoved Wattenstaafjes steriel / kleine tipp Bâtonnets ouatés stériles / petite tête Bastoncini d'ovatta sterili / testa piccola</p> <div style="display: flex; justify-content: space-between; align-items: center;"> <div style="text-align: center;"> <small>Dansu A/S Sandbakken 1-3 3660 Stenløse, Danmark</small> </div> <div style="text-align: center;"> <small>1639</small> </div> <div style="text-align: center;"> <small>Applimed SA Pra de Plan 1 CH-1618 Châtel-St-Denis</small> </div> <div style="text-align: center;"> <small>CS</small> </div> </div> <div style="text-align: center; margin-top: 10px;"> <small>(01)07640103631515(17)260327(10)210427</small> </div> </div>	
IV/3	Place of sterilisation:	ETO gas sterilization according to ISO 11135-1, By CSSR (Centre de sterilisation de Suisse Romande)	

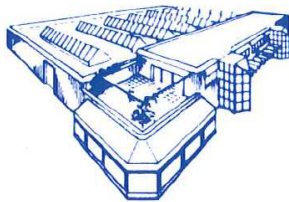


PRODUCT DATA

I-Product Description & Classification		
I/1	Designation	Wood applicator
I/2	General Description	Applicator with Wooden Stick, 150x2.5mm, Medium cotton head
I/3	GMDN N°	34954
I/4	Reference Number	1111413
I/5	Class	<input type="checkbox"/> I <input checked="" type="checkbox"/> Is <input type="checkbox"/> IIa According to Directive 93/42/EEC - Annex IX, Rules 4 and 5
I/6	Notify Body	CE ₀₁₂₀
I/7	Composition Raw Material	100 % cotton and wooden stick
I/8	Sterile	<input checked="" type="checkbox"/> Yes <input type="checkbox"/> No If, yes <input checked="" type="checkbox"/> ETO <input type="checkbox"/> Gamma <input type="checkbox"/> Steam
I/9	First Packaging	Number of pieces in Packaging: 2 pieces per peel pack
I/10	Picture	



II- Description of use and operating of device		
II/1	Intended use/indication:	A device intended used to apply a substance (eg medication) to the Ear/Nose/throat (ENT) treatment and to care the superficial wound for absorption, barrier or compressing Alone, single use and transient use (<60min)
II/2	User:	Medical staff or professional
II/3	Contraindication of use	Not use, if the device is dirty or visual defects If packaging is damaged or open do not use. Sterile unless package is damaged.
II/4	Special requirement for use	Should not be use with Alcohol
II/5	Handling and storage conditions:	Standard conditions or atmospheric conditions


III- Supplementary information		
III/1	Compatibility with others medical devices (in the set)	n/a
III/2	Device contents animal substances	NO
III/3	Device contents medicinal substances	NO
III/4	Device comes into contact with substances	Applimed has demonstrated compatibility with the following substances: - Water - Peroxide water (3%) - Betadine - Octenisept - Eosine - Merfen
III/5	Used with Accessories	NO
III/6	Used in combination with other Medical Device	NO

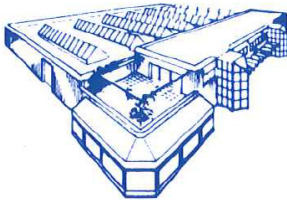

Zone industrielle
Route de Pra de Plan 1
CH-1618 CHATEL-ST-DENIS
Switzerland

TEL : +41 (0) 21 948 92 74
FAX : +41 (0) 21 948 90 36
E-mail : info@applimed.ch
Website : www.applimed.ch

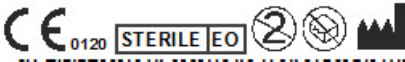IV- MANUFACTURING		
IV/1	Manufacturing Operations	<input checked="" type="checkbox"/> peel pack <input checked="" type="checkbox"/> pouches <input type="checkbox"/> boxes <input checked="" type="checkbox"/> cartons <input type="checkbox"/> cleaning <input checked="" type="checkbox"/> sterilization
		2 pieces per peel pack 100 peel packs per pouch 10 pouches per carton
IV/2	<u>Labelling: According to Directive 93/42/EEC Annex I § 13.3 and ISO 15223-1</u> <ul style="list-style-type: none"> Name and description of product Name and address of manufacturer LOT N° Article N° Year of manufacturing Expiry date Quantity Sterile Method of sterilisation Indication for single use Not use if the pack is opened CE 0120 	Example of a label Applicateur bois 150x2.5 mm REF 1111413 moyenne tête coton 2 pcs/peel pack 2018-11 100 peel packs/sachet 10 sachets/carton LOT 13.11.15   (01)07640103630433(17)181115(10)131115
IV/3	Place of sterilisation:	ETO gas sterilization according to ISO 11135-1, By CSSR (Centre de sterilisation de Suisse Romande)

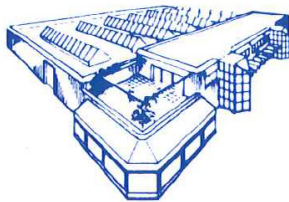


PRODUCT DATA

I-Product Description & Classification		
I/1	Designation	Wood applicator
I/2	General Description	Applicator with wooden stick, 150x2.5mm, big cotton head, continental shape
I/3	GMDN N°	34954
I/4	Reference Number	1121243
I/5	Class	<input type="checkbox"/> I <input checked="" type="checkbox"/> Is <input type="checkbox"/> IIa According to Directive 93/42/EEC - Annex IX, Rules 4 and 5
I/6	Notify Body	CE ₀₁₂₀
I/7	Composition Raw Material	100 % cotton and wooden stick
I/8	Sterile	<input checked="" type="checkbox"/> Yes <input type="checkbox"/> No If, yes <input checked="" type="checkbox"/> ETO <input type="checkbox"/> Gamma <input type="checkbox"/> Steam
I/9	First Packaging	Number of pieces in Packaging: 2 pieces per peel pack
I/10	Picture	

II- Description of use and operating of device		
II/1	Intended use/indication:	A device intended used to apply a substance (eg medication) to the Ear/Nose/throat (ENT) treatment and to care the superficial wound for absorption, barrier or compressing Alone, single use and transient use (<60min)
II/2	User:	Medical staff or professional
II/3	Contraindication of use	Not use, if the device is dirty or visual defects If packaging is damaged or open do not use. Sterile unless package is damaged.
II/4	Special requirement for use	Warning: fragile product (can break). Handle with Care Should not be use with Alcohol
II/5	Handling and storage conditions:	Standard conditions or atmospheric conditions

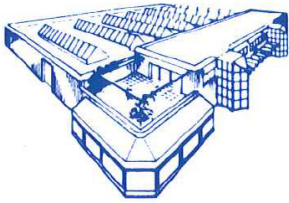

III- Supplementary information		
III/1	Compatibility with others medical devices (in the set)	n/a
III/2	Device contents animal substances	NO
III/3	Device contents medicinal substances	NO
III/4	Device comes into contact with substances	Applimed has demonstrated compatibility with the following substances: - Water - Peroxide water (3%) - Betadine - Octenisept - Eosine - Merfen
III/5	Used with Accessories	NO
III/6	Used in combination with other Medical Device	NO



Zone industrielle
 Route de Pra de Plan 1
 CH-1618 CHATEL-ST-DENIS
 Switzerland

TEL : +41 (0) 21 948 92 74
 FAX : +41 (0) 21 948 90 36
 E-mail : info@applimed.ch
 Website : www.applimed.ch

IV- MANUFACTURING

IV/1	Manufacturing Operations	<input checked="" type="checkbox"/> peel pack <input checked="" type="checkbox"/> pouches/bags <input type="checkbox"/> boxes <input checked="" type="checkbox"/> cartons <input type="checkbox"/> cleaning <input checked="" type="checkbox"/> sterilization	2 pieces per peel pack 50 peel packs per pouch 5 pouches per carton
IV/2	<u>Labelling: According to Directive 93/42 CEE Annex I § 13.3 and ISO 15223-1</u> <ul style="list-style-type: none"> Name and description of product Name and address of manufacturer LOT N° Article N° Expiry date Quantity Sterile Method of sterilisation Indication for single use Not use if the pack is opened CE 0120 	<p style="text-align: right;">Example of a label</p> <p> Applicateur bois 150x2.5 mm REF 1121243 2 Appicateurs, grosse tête continentale 50 Peelpacks/sachet 2018-11 5 Sachets/carton LOT 13.11.21 </p> <div style="text-align: center;">  <p style="font-size: small;">Applimed SA Pra de Plan 1 CH-1618 Châtel-St-Denis</p>  <p>(01)07640103633175(17)181121(10)131121</p> </div>	
IV/3	Place of sterilisation:	ETO gas sterilization according to ISO 11135-1, By CSSR (Centre de sterilisation de Suisse Romande)	